



DUAL
PLATFORM

SHOULDER
SYSTEM

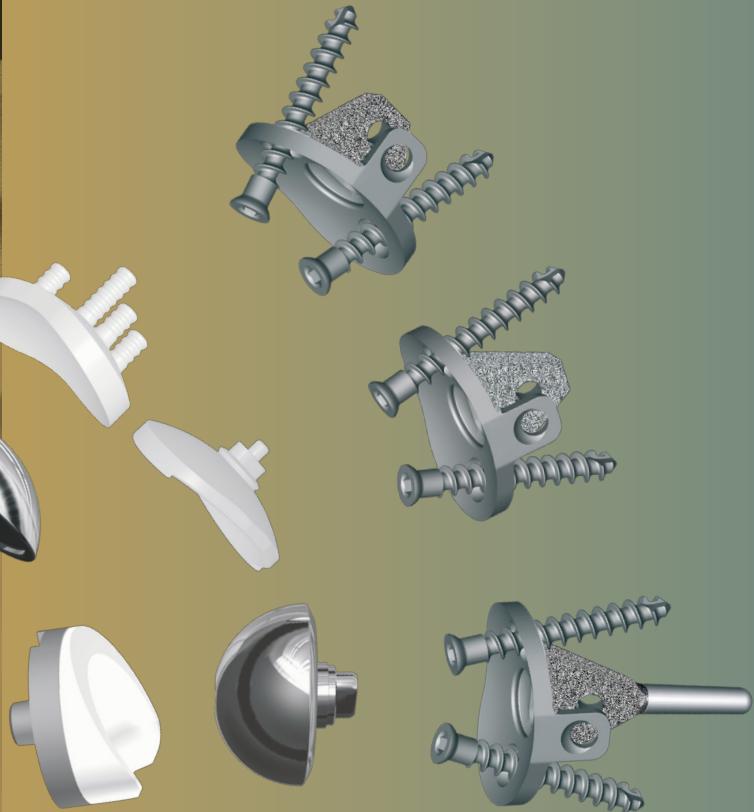

ARROW®
PRIME

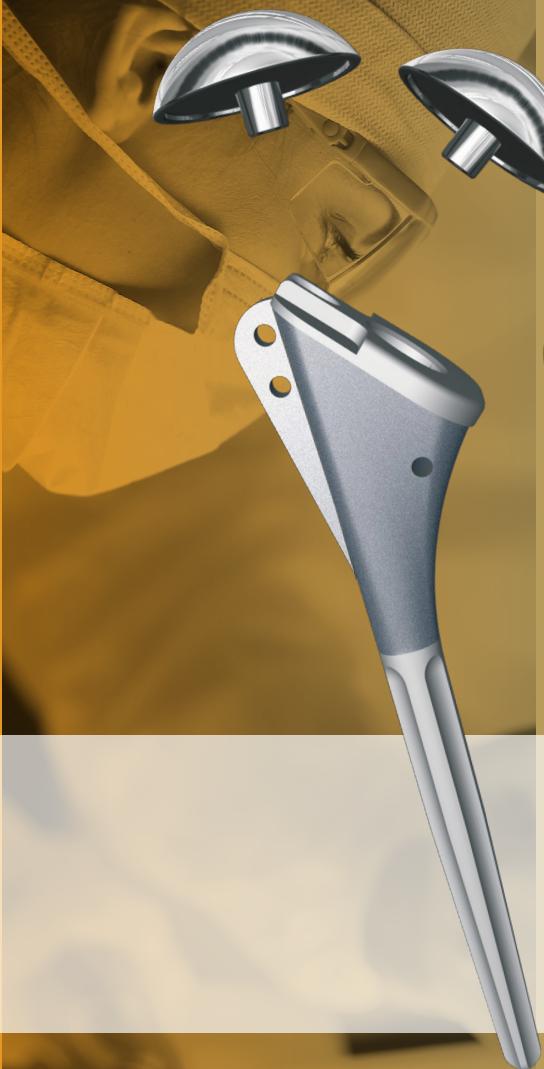


FH ORTHO™
www.FHortho.com



DUAL
PLATFORM
SYSTEM

- ARROW® PRIME has been designed to reproduce patient anatomy and reestablish soft-tissue equilibrium. Featuring a wide-range of humeral stem and head options, the ARROW® PRIME is well-suited to address various morphologies.
- The ARROW® PRIME Dual-Platform System addresses the continuum of shoulder arthroplasty indications in an efficient, four-tray instrument layout.

- 12 years of successful clinical history
- Humeral and glenoid platforms offer convertibility options
- Proven, universal metal-back glenoid baseplate
- Lateralized system optimizes motion and minimizes scapular notching

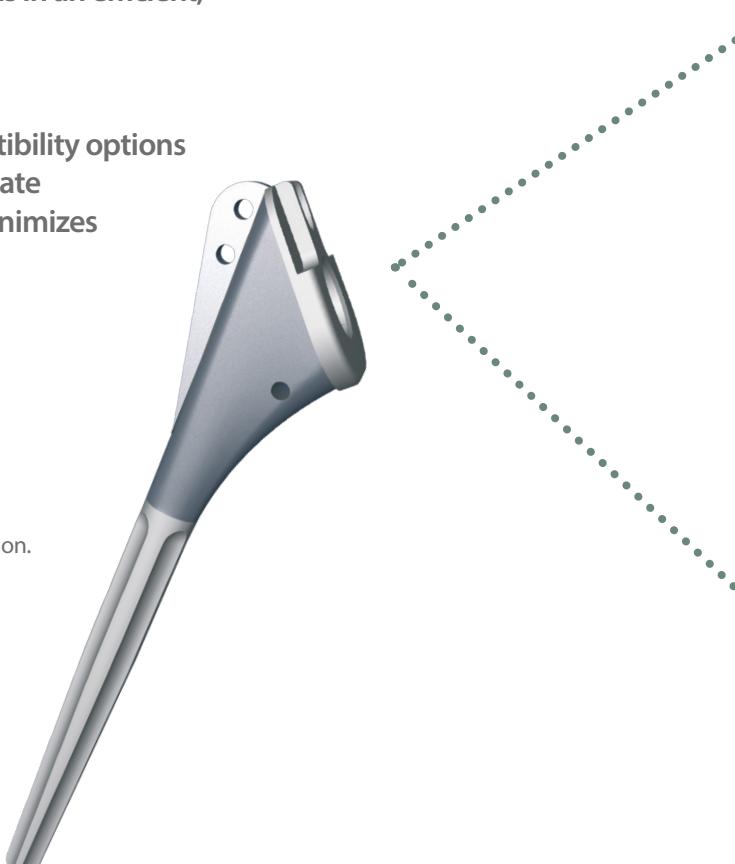

- Titanium stem features corundum-blasted metaphyseal surface to optimize fixation.
- Fixed neck-angle of 135 degrees.
- Stability provided by metaphyseal volume and lateral fin.
- Humeral stem can be implanted with or without cement.
- Stem diameters of 6, 8, 10, 12, 14 and 16 mm.
- Proportional stem length limits varus or valgus malpositioning.
- Female Morse-Taper facilitates glenoid exposure and visualization.
- The lateral fin suture holes support tuberosity reduction and fixation for fracture reconstruction.



ECCENTRIC HEAD



CONCENTRIC HEAD

ANATOMIC



CEMENTED GLENOID

- Pegged, curved-back polyethylene implant conserves bone and maximizes stability.
- A 4mm "mismatch" allows for proper humeral head translation and reduces the risk of glenoid loosening.




GLENOID INSERT

- Consistent 4mm thick polyethylene insert reduces contact stress and improves wear properties.
- Secure, single-step locking of insert to porous glenoid baseplate.

HUMERAL HEADS

- CoCr humeral head.
- Concentric and eccentric heads are available to optimize medial and posterior offsets.
- Diameters of 40 to 54mm available to achieve ideal humeral coverage.
- Multiple head-heights available to restore soft-tissue tensioning.

REVERSE

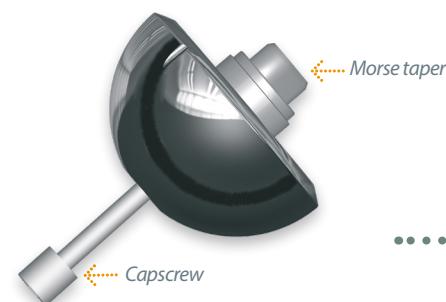


HUMERAL INSERT

- Patented medial cut-out minimizes scapular notching.
- Increased congruity to optimize stability.
- Pre-assembled humeral insert prevents dissociation.

GLENOSPHERE

- Three diameters available (36, 39, and 42mm).
- Morse taper couples glenosphere to baseplate.
- Capscrew secures glenosphere to the baseplate, preventing dissociation.

HEMI

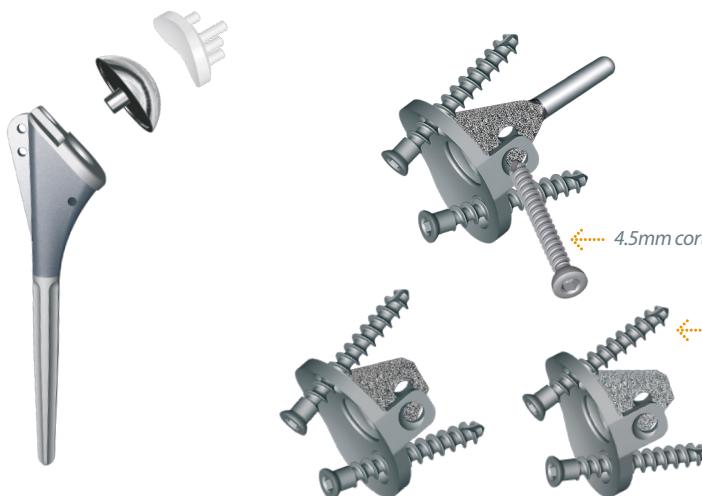


4.5mm cortical screw (optional)*

5.5mm cancellous screw



CEMENTLESS POROUS GLENOID BASE

- Titanium porous coating.
- Multiple sizes available to accommodate variable anatomy.

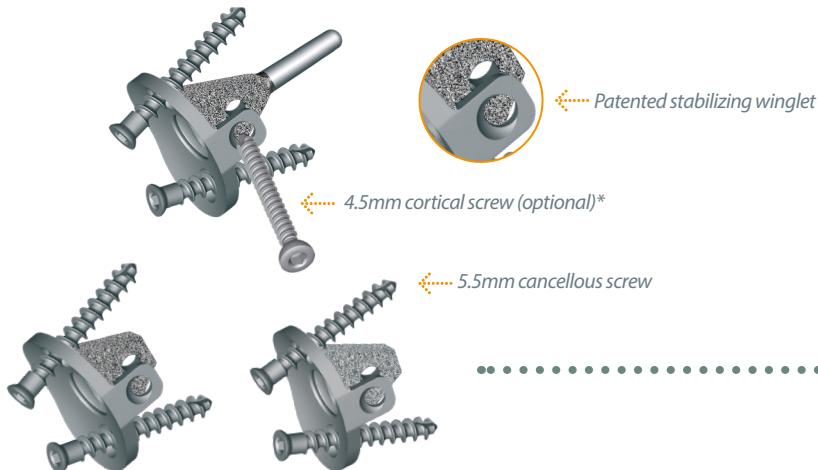

Patented stabilizing winglet

4.5mm cortical screw (optional)*

5.5mm cancellous screw



CEMENTLESS POROUS GLENOID BASE

- Divergent 5.5mm cancellous bone screws.
- Multiple sizes available to accommodate variable anatomy.
- Curved-back design improves stability and force transfer.

* OPTIONAL
Cortical bone screw for compromised glenoid bone and for fractures of the glenoid

REFERENCE NUMBERS

HUMERAL STEM



REFERENCE	DIAMETER	HEIGHT
267 360	Ø 06 (CEMENTED ONLY)	100
265 102	Ø 08	120
265 103	Ø 08	170
265 104	Ø 10	125
265 105	Ø 12	130
265 106	Ø 14	135
267 361	Ø 16	140

CONCENTRIC HUMERAL HEAD



REFERENCE	DIAMETER	HEIGHT
265 107	Ø 40	15
265 108	Ø 40	17
265 109	Ø 44	16
265 110	Ø 44	18
265 111	Ø 46	16
265 112	Ø 46	18
265 113	Ø 46	21
265 114	Ø 48	16
265 115	Ø 48	18
265 116	Ø 48	21
265 117	Ø 50	17
265 118	Ø 50	19
265 119	Ø 50	21
265 120	Ø 52	19
265 121	Ø 52	21
265 122	Ø 54	19
265 123	Ø 54	21

ECCENTRIC HUMERAL HEAD



REFERENCE	DIAMETER	HEIGHT
265 124	Ø 44	16
265 125	Ø 44	18
265 126	Ø 46	16
265 127	Ø 46	18
265 128	Ø 46	21
265 129	Ø 48	16
265 130	Ø 48	18
265 131	Ø 48	21
265 132	Ø 50	17
265 133	Ø 50	19
265 134	Ø 50	21
265 135	Ø 52	19
265 136	Ø 52	21

CEMENTED GLENOID



REFERENCE	SIZE
265 137	44
265 138	46
265 139	48
265 140	50

POROUS GLENOID IMPLANT



REFERENCE	SIZE
267 702	44S
267 701	44
267 704	46
267 705	48
268 698	44S-LP*
267 703	44-LP*
268 699	46-LP*

* LONG POST

GLENOID INSERT



REFERENCE	SIZE
265 157	44
265 158	46
265 159	48

GLENOSPHERE



REFERENCE	DIAMETER
265 150	Ø 36
265 151	Ø 39
265 152	Ø 42

HUMERAL INSERT



REFERENCE	DIAMETER	HEIGHT
265 141	Ø 36	00
265 142	Ø 36	05
265 143	Ø 36	10
265 144	Ø 39	00
265 145	Ø 39	05
265 146	Ø 39	10
265 147	Ø 42	00
265 148	Ø 42	05
265 149	Ø 42	10

→ special sizes for larger anatomies, trauma or revision cases

CANCELLOUS BONE SCREW - sterile -



REFERENCE	DIAMETER	LENGTH
265161	Ø 5.5	24
265162	Ø 5.5	28
265163	Ø 5.5	32
265164	Ø 5.5	36
265165	Ø 5.5	40
265166	Ø 5.5	45
265167	Ø 5.5	50

CORTICAL BONE SCREW - sterile -



REFERENCE	DIAMETER	LENGTH
265 168	Ø 4.5	32
265 169	Ø 4.5	34
265 170	Ø 4.5	36
265 171	Ø 4.5	38
265 172	Ø 4.5	40

BIBLIOGRAPHY

- The main cause of instability after unconstrained shoulder prosthesis is soft tissue deficiency
Jean Kany MD, Jijo Jose DNB, Denis Katz MD, Jean David Wertheil MD, Padmanaban Sekaran, MSc PT, Rajkumar S. Amaravathi DNB, Philippe Valenti MD - J Shoulder Elbow Surg (2017)
- Regarding Metal-backed glenoid implant with polyethylene insert is not a viable long-term therapeutic option
Denis Katz - J Shoulder Elbow Surg (2016)
- Does lateralisation of the centre of rotation in reverse shoulder arthroplasty avoid scapular notching? Clinical and radiological review of one hundred and forty cases with forty five months of follow-up (2015)
Denis Katz & Philippe Valenti & Jean Kany & Kamil Elkholi & Jean-David Wertheil - International Orthopaedics (SICOT) - DOI 10.1007/s00264-015-2976-3
- Uncemented metal-back glenoid component in revision of aseptic glenoid loosening: a prospective study of 10 cases with a minimum follow-up of 2 years.(2014)
Valenti P, Valle P, Sauzieres P, Boughebri O, Moraiti T, Maquès A, Spies CK, Unglaub F, Lecâtre FM. Arch Orthop Trauma Surg; 2014 Jun;134(6):785-91.
- A convertible shoulder system : is it useful in total shoulder arthroplasty revisions (2014)
J. Kany, T. Amouyal, O. Flaman, D.Katz, P. Valenti : International orthopaedics DOI 10.1007/s00264-014-2563-z
- New design of a cementless glenoid component in unconstrained shoulder arthroplasty : a prospective medium-term analysis of 143 cases (2013)
Eur J Orthop Surg Traumatol (2013). 23 :27-34 - DOI 10.1007/s00590-012-1109-6
- Results of 30 reverse shoulder prostheses for revision of failed hemi- or total shoulder arthroplasty.(2013)
Philippe Valenti, Alexandre Sahin Kilinc, Philippe Sauzieres, Denis Katz European Journal of Orthopaedic Surgery & Traumatology (Impact Factor: 0.1). 10/2013; DOI:10.1007/s00590-013-1332-9
- How to deal with glenoid type B2 or C? How to prevent mistakes in implantation of glenoid component? (2012)
Jean Kany, Denis Katz - Eur J Orthop Surg Traumatol. 2013 May;23(4):379-85. Epub 2012 Nov 2.
- Mid-term outcome of reverse shoulder prosthesis in complex proximal humeral fractures (2012)
P. Valenti , D.Katz, A.Kilinc , K.Elkholti , V.Gasiunas Acta Orthop Belg (2012) 78-4
- Lateralization in RSA-metallic or bony lateralization ?(2012)
P. Valenti , D.Katz, P.Sauzieres , J.Kany , A.Kilinc Shoulder arthroscopy and arthroplasty –Current Concepts 2012 tome 2
- The case for the metal-back glenoid design in TSA (2012)
D.Katz , P.Sauzieres,P Valenti J.Kany Eur J Orthop Surg Traumatol(2012) 22 :9-16
- Do Less Medialized Reverse Shoulder Prostheses Increase Motion and Reduce Notching? (2011)
Valenti P, Sauzieres P, Katz D, Kalouche I, Kilinc AS. Clin Orthop Relat Res. (2011) Mar 15. - DOI 10.1007/s11999-011-1844-8
- Reverse shoulder arthroplasty: does reduced medialisation improve radiological and clinical results?(2009)
Kalouche I, Sevivas N, Wahegaonker A, Sauzieres P, Katz D, Valenti P - Acta Orthop Belg. 2009 Apr;75(2):158-66.
- The reverse shoulder prosthesis surgical technique.(2008)
Valenti P, Sauzieres P, Cogswell L, O'Toole G, Katz D. Tech Hand Up Extrem Surg. 2008 Mar;12(1):46-55.
- A history of the reverse shoulder prosthesis (2007)
D.Katz, O'Toole G, Cogswell L,P.Valenti ,P.Sauzieres - Int J Shoulder Surg 2007, 1 :108-13
- 20 Cases of chronic anterior and posterior dislocations of the shoulder treated by arthroplasty. Considerations about prognosis and treatment options.(2007)
D.Katz , P.Valenti,P. Sauzieres - Poster 10 th International Congress of Shoulder Surgery(ICSS)sept 17-20 2007-Bahia-BRASIL
- 20 cases of chronic dislocations of the shoulder treated by arthroplasty.Pronostic and therapeutic implications (fr)(2007)
D.Katz,P Valenti ,A El Hadi - Rev Chir Orthop Oct 2004,Vol 90,Sup 6, pp 25145