


# FH LOOP

NON-ABSORBABLE FIBER LOOP SUTURE

**FH LOOP is a non-absorbable fiber loop suture composed of a braided ultra-high molecular weight polyethylene (UHMWPE/Polyester) loop coated with uniform polybutylene adipate.**

- Thread:  $\varnothing$ 0.50 (USP 2)
- Length: 60 cm (loop length : 50 cm)
- Colors: blue and white
- 300-series steel (steel and inox alloy) straight, silicone coated and Cobra Black® (black) needle, 60 mm length

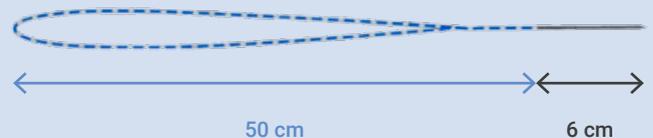

## Product characteristics

FH LOOP needle is composed of 300-Series stainless steel to enhance their penetration force, strength rate and provide good ductility avoiding risks of breakage.

FH LOOP fiber loop suture integrates an inter-locking core technology which provides all fiber configurations with a sturdy core within center of the structure to act as backbone for better knot tying and handling characteristics.

FH LOOP suture elicits a minimal inflammatory reaction in tissues, followed by gradual encapsulation of the suture by fibrous connective tissue.

The FH LOOP suture is non-absorbable and there are no significant changes in tensile strength retention known to occur in vivo.

## Indications

FH LOOP is indicated for use in the fixation of bone and soft tissue in orthopedic procedures requiring ligament or tendon reconstruction.

This fiber loop suture are highly appropriate for tendinous graft preparation, ACL repair, rotator cuff repair, Achilles tendon repair, other soft tissue and anchor applications.

## Packaging

FH LOOP is supplied sterile, sterilized by ethylene oxide including a black needle for less glare and higher visibility.


Do not re-sterilize.

# FH LOOP

NON-ABSORBABLE FIBER LOOP SUTURE



**E-braid technology**  
More fibers for better strength

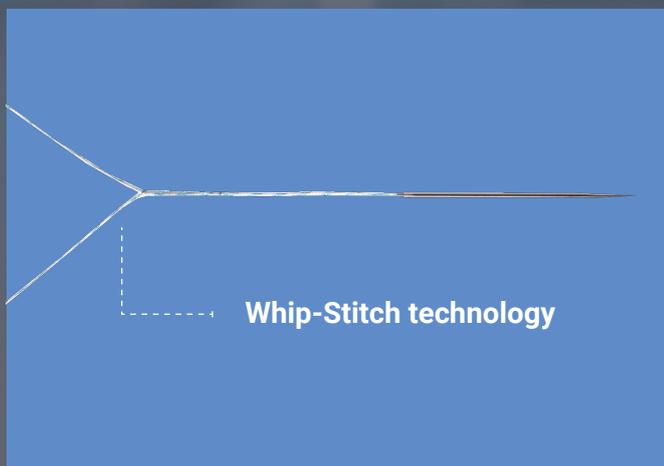

**Whip-Stitch technology**

## Reference

FH ref.	Supplier ref.	Designation
271 197	HS173	FH LOOP fiber loop suture

## GROUP FH ORTHO™

### FR, FH ORTHO SAS

3 rue de la Forêt  
68990 HEIMSBRUNN - FRANCE  
Tél. +33 (0)3 89 81 90 92 / Fax : +33 (0)3 89 81 80 11  
info@fhortho.com / www.fhortho.com

### USA, FH ORTHO INC.

4908 N. Elston, Chicago, Illinois 60630 - USA  
Tel.: +1 (773) 290 1039 / 844-77 FHINC  
Fax: +1 (708) 667 7618  
info-us@fhorthopedics.com / www.fhortho.com

### UK, FH ORTHO LTD

Unit 1b, Century Park, Valley way  
Swansea Enterprise Park, Swansea, SA6 8RP - UK  
Tel.: +44 (0) 1792 464792  
Fax: +44 (0) 844 412 7674  
customer-servicesUK@fhorthopedics.com  
www.fhortho.com

### PL, FH ORTHO POLSKA

Ul. Garbary 95/A6,  
61-757 Poznan - POLSKA  
Tel: +48 61 863 81 27 / Fax: +48 61 863 81 28  
biuro@implants24.pl  
www.fhortho.com